

# **BioTek Synergy H1 Microplate Reader for Absorbance Measurements**

Teng-Jui Lin

Department of Chemical Engineering, University of Washington

**Surface and Colloid Science**

# Beer-Lambert law relates absorbance and solution concentration linearly


- **Intensity**  $I$  - measure of how much light
- **Transmittance**  $T$  - ratio of transmitted light vs. incident light at a wavelength  $\lambda$

$$T = \frac{I}{I_0} \quad \frac{\text{final}}{\text{initial}} \leq 1$$

- **Absorbance**  $A$  - capacity of a substance to absorb light at a wavelength  $\lambda$

$$A = -\log_{10}(T) = -\log_{10}\left(\frac{I}{I_0}\right) = \log_{10}\left(\frac{I_0}{I}\right) \in (0, 2)$$



$\geq 1$   
magnitude changes

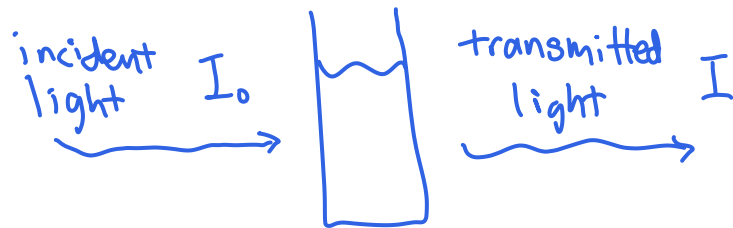

- **Beer-Lambert law** - absorbance varies linearly with solution concentration and path length



$$A = \epsilon b C$$

- $\epsilon$  - extinction coefficient  $\text{const}$  [
- $b$  - path length [cm]
- $C$  - solution concentration [M, mM]

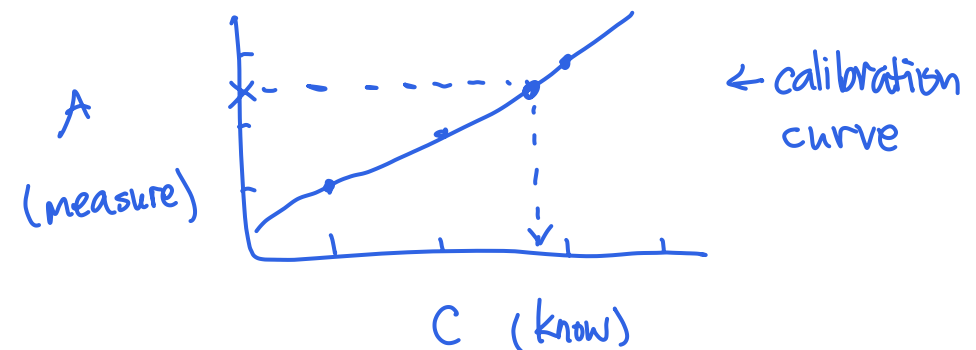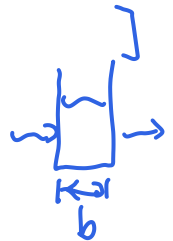


# 96-well plates contains 200 $\mu\text{L}$ sample of interest in each well

- 200  $\mu\text{L}$  sample in each well
- At least one blank control
- Avoid air bubble (reverse pipetting)
- Label the wells

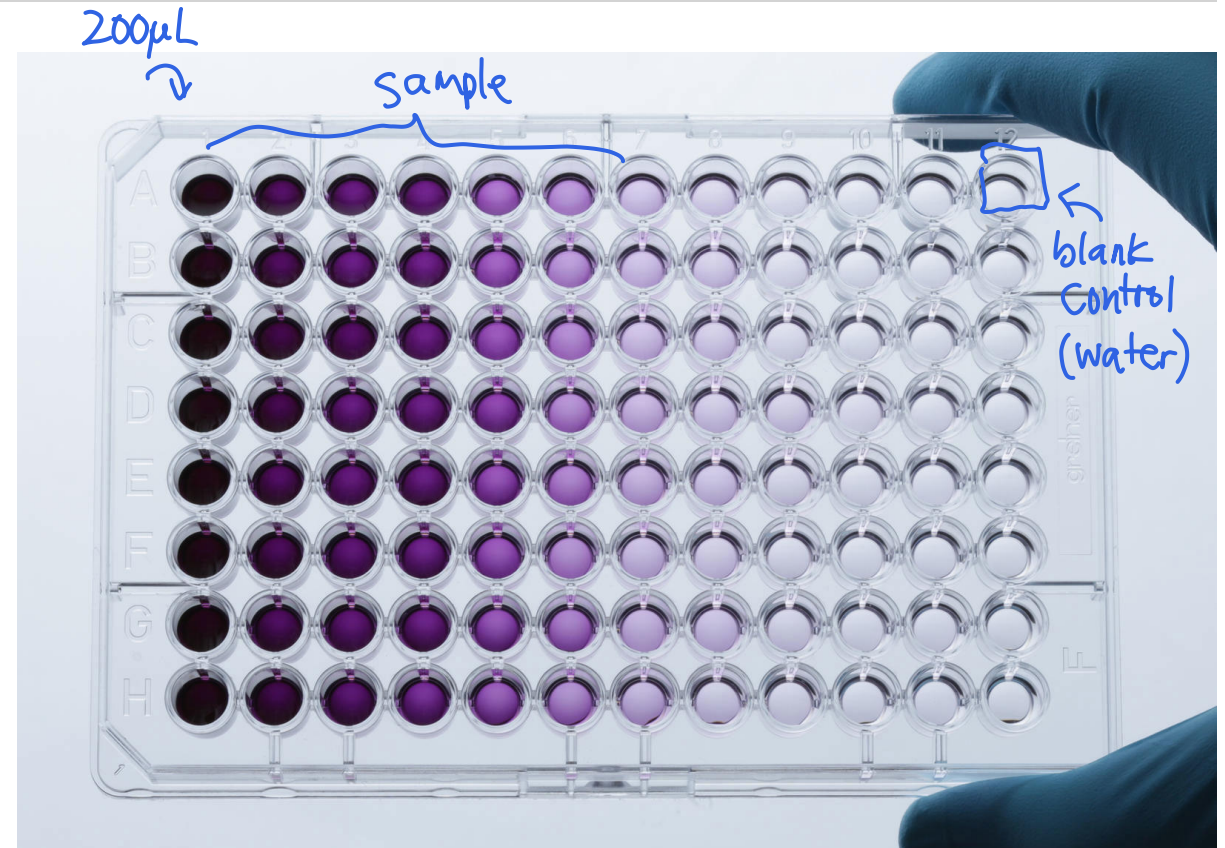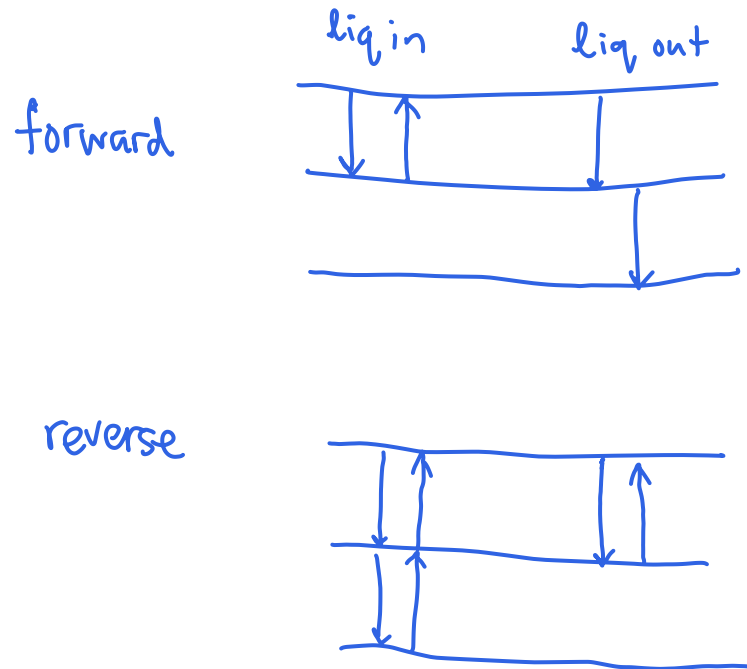

# Microplate reader location, access, and training

- The microplate reader (BioTek Synergy H1) is located in the Bindra Innovation Laboratory (Benson Hall 121).
- Book usage time on shared Google Calendar
- Log usage time on logbook

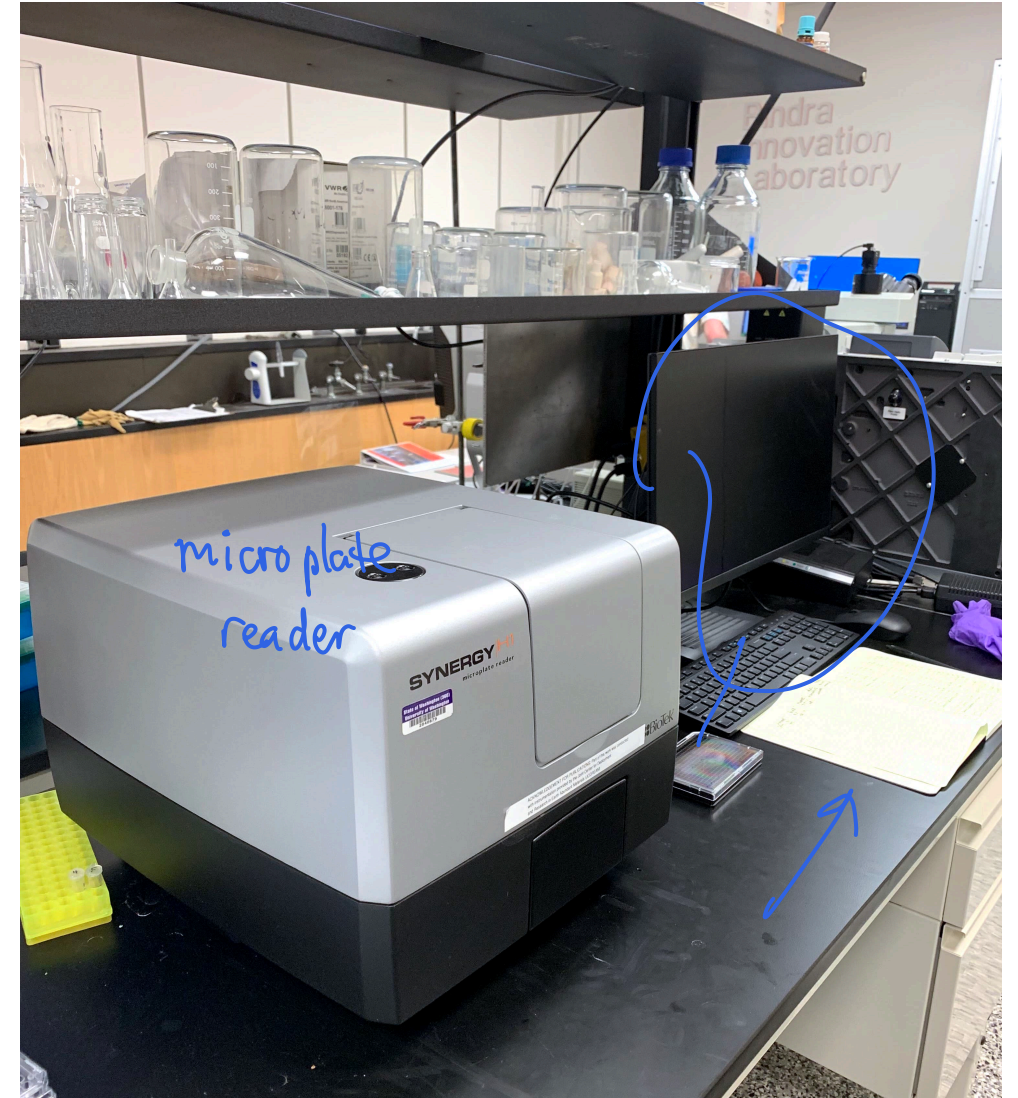



# Microplate reader startup

- Turn on the microplate reader
  - Wait for self-diagnosis
  - Create empty Google Sheet
- Open "Gen5 3.09" software




# Plate reading settings - absorbance wavelength

---

- Measurement mode (default)
  - Read method: Absorbance
  - Read type: Endpoint/Kinetic
  - Optics type: Monochromators
- Absorbance wavelength
  - Lab 3-1: CMC by Dye Titration
    - $\lambda = 615$  nm for pinacyanol chloride
  - Lab 3-3: Aggregation of Clay
    - $\lambda = 860$  nm for turbidity
  - Lab 4-4: Bubble Fractionation
    - $\lambda = 590$  nm for crystal violet

# Plate reading settings - plate layout

- Plate Layout
  - Select “Blanks” and “Samples”

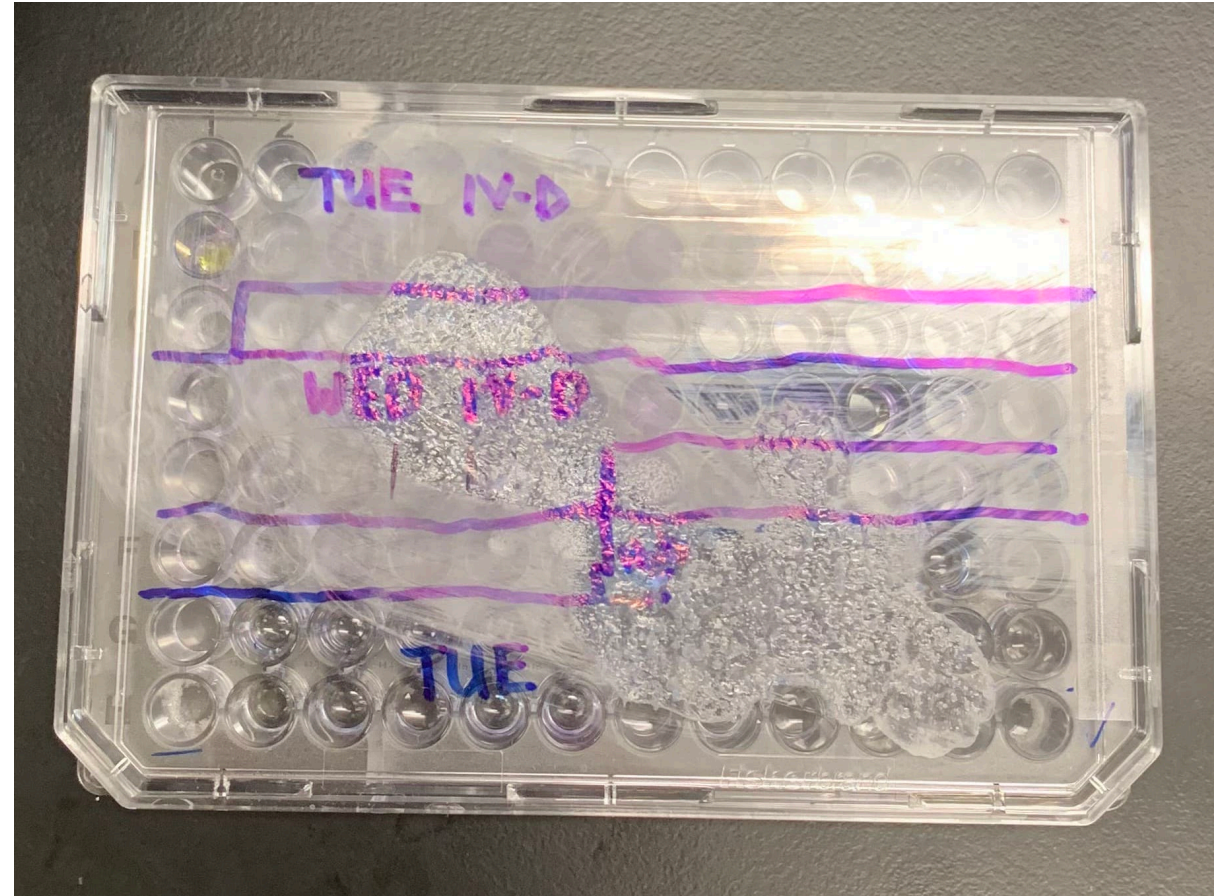



# Absorbance measurement

- Place the plate into the plate holder
  - A1 well is on the top right.
- Read the plate using computer software
- Remove the plate from the plate holder.





## Data recording

---

- Export the data in both matrix and stats form
  - Use blank-subtracted absorbance

# Microplate reader shutdown

---

- Close all programs on the computer.
- Push the IN/OUT button for the plate holder so the plate holder retracts.
- Push the ON/OFF button for the plate reader so the instrument is turned off.

# Absorbance measurements for Surface and Colloid Science Laboratory

- Lab 3-1: CMC by Dye Titration
  - $\lambda = 615$  nm for pinacyanol chloride
  - Look for changes in linear trend of absorbance
- Lab 3-3: Aggregation of Clay
  - $\lambda = 860$  nm for turbidity
  - Need calibration curve
- Lab 4-4: Bubble Fractionation
  - $\lambda = 590$  nm for crystal violet
  - Need calibration curve

