

BET (Brunauer-Emmett-Teller) Method for Surface Area Determination

Teng-Jui Lin

Department of Chemical Engineering, University of Washington

Surface and Colloid Science

BET isotherms empirically fits adsorption isotherm data

- BET method - procedure for determining surface area of finely divided solid by measurements of low-temperature gas adsorption
- BET isotherm

$$\frac{V}{V_m} = \frac{C(p/p_0)}{(1 - p/p_0)[1 + (C - 1)(p/p_0)]}$$



- Linear form of BET isotherm

$$\frac{(p/p_0)}{V(1 - p/p_0)} = \frac{1}{CV_m} + \frac{(C - 1)(p/p_0)}{CV_m}$$

- Single point method: $C \gg 1$, so intercept ≈ 0

$$\text{Volume: } V_m = V(1 - p/p_0)$$

$$\text{Mass: } X_m = X(1 - p/p_0)$$

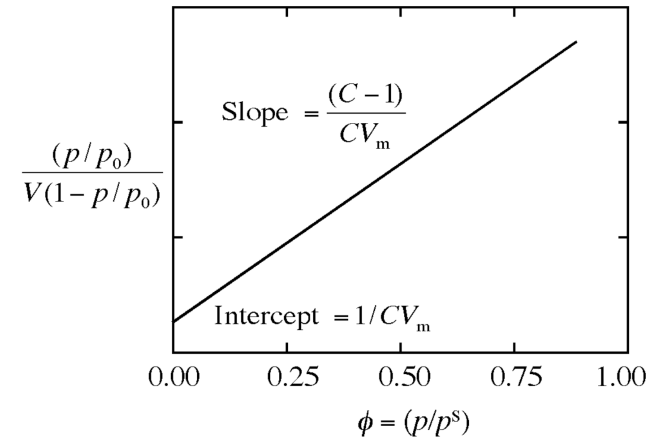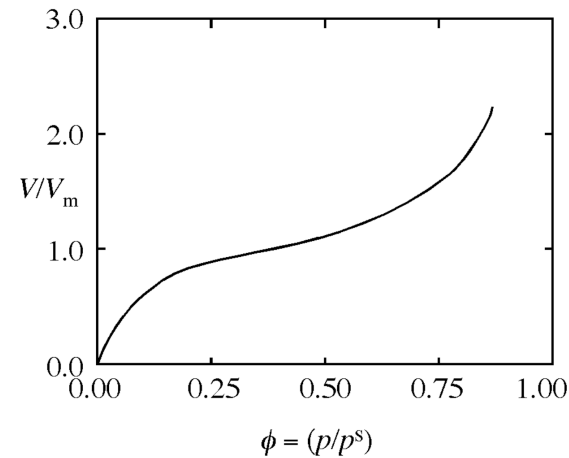

BET method allows surface area determination

- Volume of adsorbate desorbed by sample

$$V_{\text{des}} = V_c \frac{A}{A_c}$$

- Mass of adsorbate desorbed

$$X = \frac{p_a V_{\text{des}}}{RT} M_a$$

- Mass of adsorbate desorbed when solid is covered by a single adsorbate monolayer

$$X_m = \frac{p_a V_c (A/A_c)}{RT} M_a (1 - p/p_0)$$

- Total surface area

$$S_t = X_m \frac{N_0}{M_a} A_{cs}$$

- Specific surface area

$$\boxed{\Sigma = \frac{S_t}{m}}$$

Initial configuration of FlowSorb Surface Area Analyzer

- Initial configuration
 - Power: OFF
 - Toggle valve: open
 - Sample holders are installed at DEGAS, TEST, and COLD TRAP locations
 - Insulating container is installed at sample holder for COLD TRAP



Startup procedure of FlowSorb Surface Area Analyzer

- Startup
 - Open inert gas flow (30% nitrogen, 70% helium), wait 5 min
 - Adjust flow meter to the calibration mark
 - Pour liquid nitrogen into an insulating container, and pour the liquid nitrogen to the insulating container at COLD TRAP, wait for 10 min
 - Power switch: ON, wait for 30 min

Preparing activated carbon sample for BET measurement

- Sample preparation
 - Weigh 0.02 g activated carbon in the sample holder
 - Place the sample holder on DEGAS and heat with a heating mantle, wait for 15 min (concurrent with start up)

Calibrating FlowSorb Surface Area Analyzer

- Calibration
 - Flush gas syringe with evaporated nitrogen gas above the liquid nitrogen at DEGAS
 - Fill the gas syringe with 1 mL nitrogen gas
 - Wipe needle tip free of frost, set aside needle to equilibrate to room temperature
 - Zero the instrument display with course and fine zeros. Switch to surface area (SA), and clear SA display
 - Insert needle at INJECT, inject nitrogen gas at moderate rate, and withdraw
 - Wait until reading is stabilized, then calibrate the instrument to **SA = 2.84**
 - Confirm calibration by repeated injection, if necessary

Measuring surface area of activated carbon with the BET method (adsorption)

- Adsorption measurement
 - Exchange sample holders between DEGAS and TEST, so the sample is at TEST
 - Wait until readings stabilize and clear SA display
 - Pour liquid nitrogen into an insulating container
 - Place the insulating container at TEST by immersing the sample holder with liquid nitrogen. Secure the container by flipping on the container holder
 - Nitrogen gas starts to adsorb to activated carbon
 - Wait until reading is stabilized, **record the value of SA** = adsorption surface area

Measuring surface area of activated carbon with the BET method (desorption)

- Desorption measurement
 - Remove the insulating container by putting down the container holder
 - Nitrogen gas starts to desorb from activated carbon
 - Wait until reading is stabilized, **record the value of SA** = desorption surface area

Shutdown procedure of FlowSorb Surface Area Analyzer

- Shutdown
 - Power: OFF
 - Remove the activated carbon from the sample holder, and place the sample holder back on TEST
 - Turn off the inert gas flow

Calculating the specific surface area of activated carbon

- Measured

S_t : Adsorption SA, desorption SA

m : Mass of activated carbon

- Calculated

- $\Sigma = \frac{S_t}{m}$ - specific SA for adsorption and desorption

Usually use desorption data to report Σ