

# **Surface Tension Measurements by Detachment and Partial Immersion Methods**

Teng-Jui Lin

Department of Chemical Engineering, University of Washington

**Surface and Colloid Science**

# Surface tension can be measured from force methods

- Detachment method (du Noüy ring)
  - Force required to pull a solid completely through a fluid interface
- Partial immersion method (Wilhelmy slide)
  - Force required to maintain the position of a solid which penetrates a fluid interface
- Force balance
  - Down = Weight - Buoyancy + Surf. Tension

$$F_{\downarrow} = F_g - F_b + F_{\sigma}$$


*perimeter*  
↓

$$F_{\downarrow} = mg - \rho g V_{\text{disp}} + P\sigma$$


Assumes: uniform  $\sigma$ , fully wetted  $\theta = 0$

du Noüy ring



Overall view

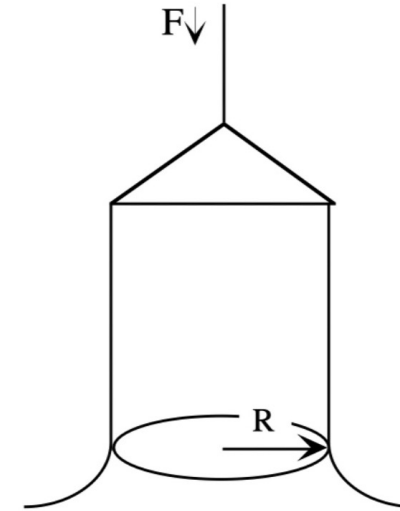

Side view

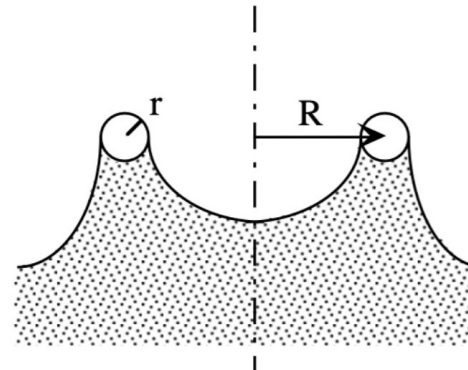

Wilhelmy slide

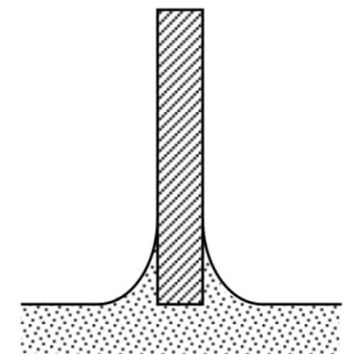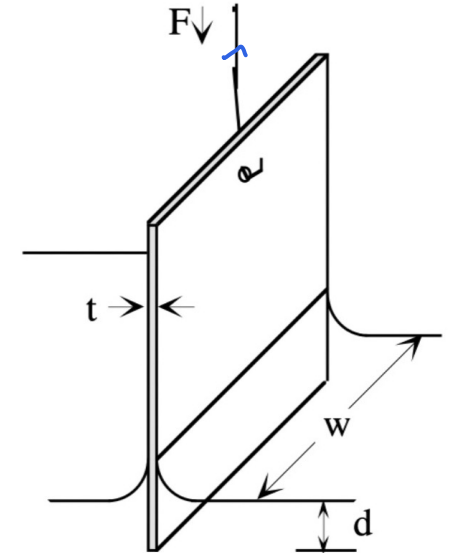

# du Noüy ring can be used to measure $\sigma$ by the detachment method

- Force balance

$$\begin{aligned}
 F_{\downarrow} &= mg - \rho g V_{\text{disp}} + P\sigma \\
 &= mg - 0 + 2\pi[(R+r) + (R-r)]\sigma \\
 &= mg + 4\pi R\sigma
 \end{aligned}$$

- Uncorrected surface tension  $\sigma^*$

$$\sigma^* = \frac{F_{\downarrow} - mg}{4\pi R}$$



- Correction factor  $F$  [cgs unit]

$$F_{\downarrow} = mg + \frac{4\pi R\sigma}{F}$$

$$F = 0.725 + \sqrt{\frac{0.01425\sigma^*}{(2\pi R)^2(\rho_l - \rho_{\text{air}})}} + 0.04534 - \frac{1.679}{R/r}$$

$$\sigma = \frac{F_{\downarrow} - mg}{4\pi R} F$$

du Noüy ring

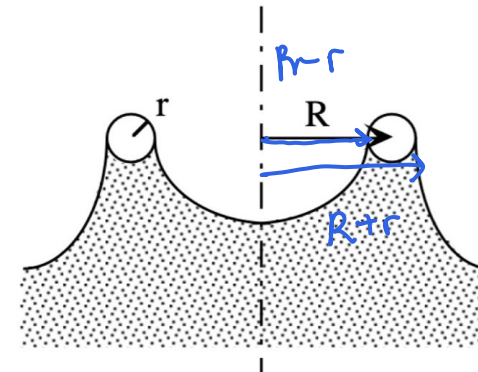
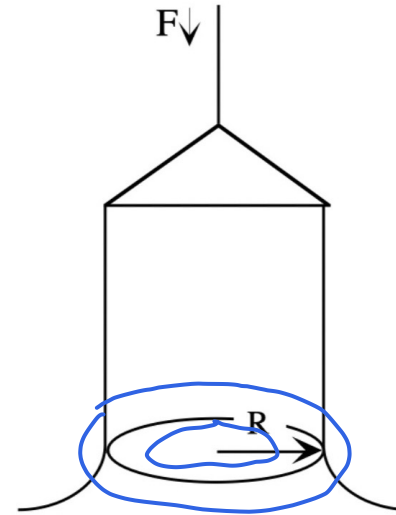


# Wilhelmy slide can be used to measure $\sigma$ by the partial immersion method


- Force balance ( $d = 0, t \ll w$ )

$$\begin{aligned} F_{\downarrow} &= mg - \rho g V_{\text{disp}} + P\sigma \\ &= \cancel{m}g - \cancel{\rho}g\cancel{t}w\cancel{d} + 2(w + \cancel{t})\sigma \\ &= mg + 2w\sigma \end{aligned}$$

Wilhelmy slide

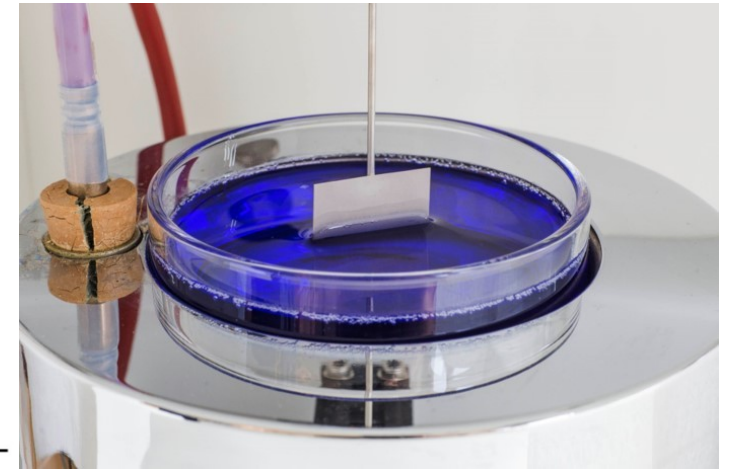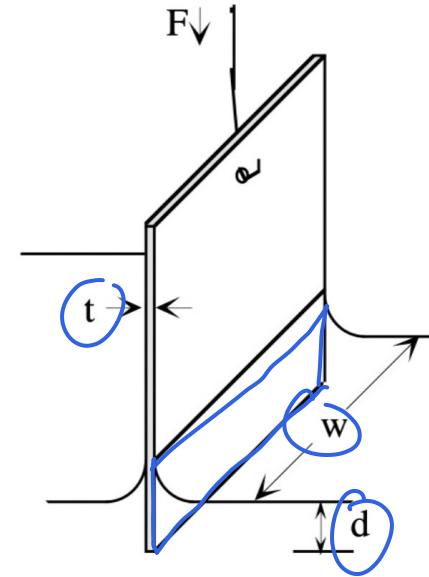

$$V_{\text{disp}} = twd$$

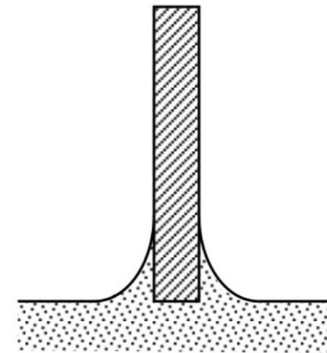

# Tensiometer is used to measure force on du Noüy rings and Wilhelmy slides

$$F_{\downarrow} \Rightarrow \sigma$$

